

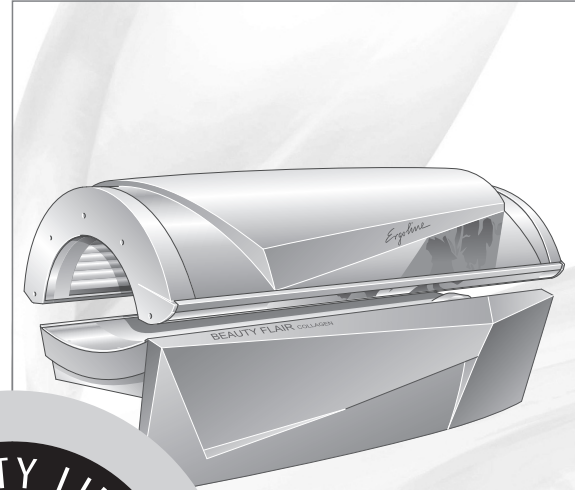
BEAUTY FLAIR

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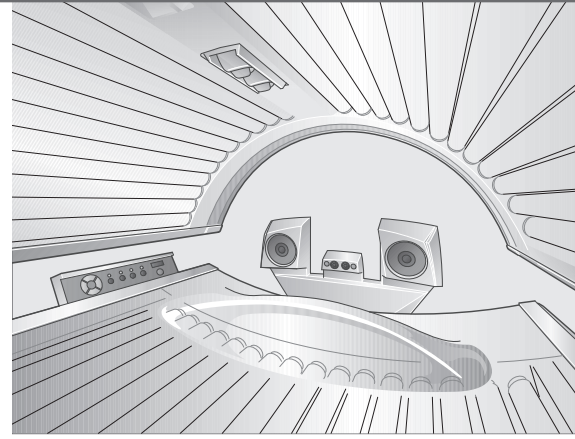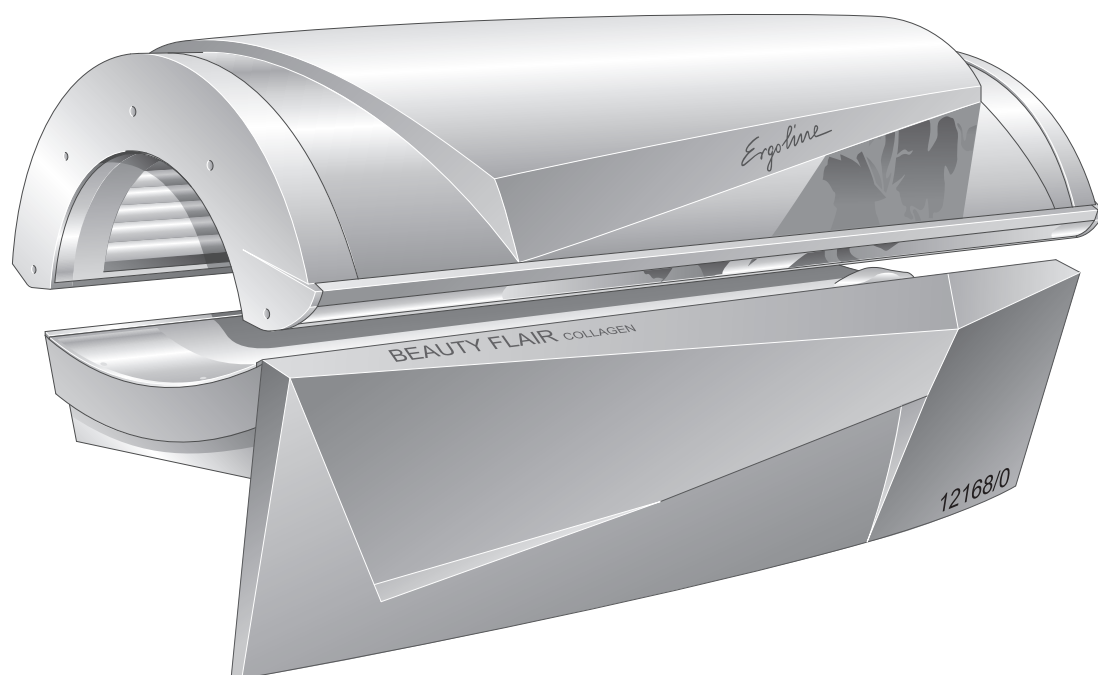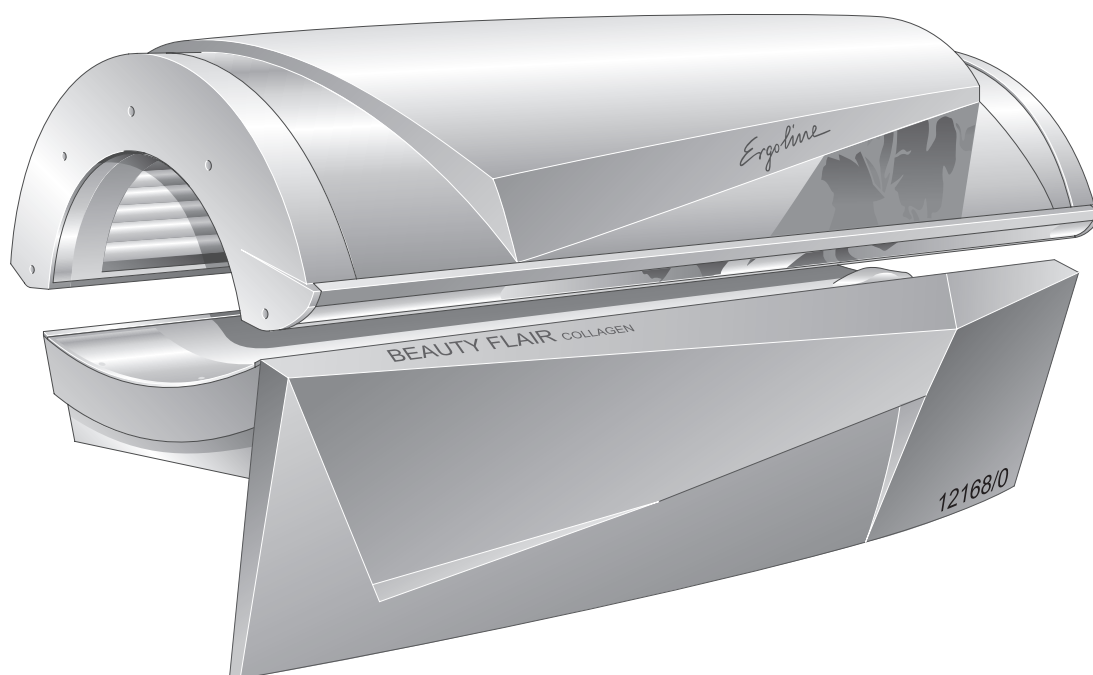
Operating instructions

Translation of the original operating instructions

BEAUTY FLAIR COLLAGEN

1008571-00 / en / 10.2009

Ergoline

Imprint

Manufacturer: JK-Products GmbH

Köhlershohner Straße
53578 Windhagen
GERMANY

Tel.: +49 (0) 22 24 / 818-140

Fax: +49 (0) 22 24 / 818-166

**Customer service /
Technical service (Spare
part orders for
components): JK-Global Service GmbH**

Rottbitzer Straße 69
53604 Bad Honnef (Rottbitze)
GERMANY

Tel.: +49 (0) 22 24 / 818-861

Fax: +49 (0) 22 24 / 818-205

E-Mail: service@jk-globalservice.de

**Spare part orders for
consumables: JK-Licht GmbH**

Eduard-Rhein-Str. 3
53639 Königswinter
GERMANY

Tel.: +49 (0) 22 44 / 9239-0

Fax: +49 (0) 22 44 / 9239-24

E-Mail: info@jk-licht.de

DANGER!

Failure to observe these instructions:

- can result in serious injury and death,
- can result in damage to the device and the environment.
- Read these instructions carefully before putting this collagen device into operation.
- Please observe the instructions and the code of conduct required to operate the device safely.
- Please make the instructions and additional information from the manufacturer available to the personnel at their workplace.

Copyright

The copyright remains with JK-Holding GmbH.

The contents may not be partially or completely duplicated, distributed or utilised for competition purposes without authorisation or disclosed to other parties.

The right to make technical changes in regard to the figures and information contained in this instruction manual is reserved!

Table of contents

1	General.....	5
1.1	Definitions	5
1.2	Symbol explanation.....	6
2	Information for the operator.....	7
2.1	Intended use	7
2.2	Foreseeable misuse	7
2.3	Operator's obligations	7
2.4	Personnel training	8
2.5	Scope of delivery.....	8
2.6	Transport, assembly and setup.....	9
2.6.1	Assembly	9
2.6.2	Ventilation	9
2.6.3	Setup location	10
2.7	Commissioning.....	10
2.8	Operation and maintenance.....	10
2.9	Decommissioning.....	11
2.10	Storage	11
2.11	Disposal	11
2.11.1	Declaration – The JK Corporate Group.....	11
2.11.2	Environmental regulations – disposal of lamps and batteries	11
2.11.3	Packaging	12
2.11.4	Disposal of old devices	12
2.12	Directives	12
2.13	Export.....	12
2.14	MP3 music (optional)	13
2.15	Instruction signs on the device	13
2.16	Technical modifications.....	13
2.17	Manufacturer's guarantee	14
2.18	Warranty and liability exclusions	14
3	Description	15
3.1	Scope of delivery.....	15
3.1.1	Optional.....	15
3.2	Equipment.....	15
3.3	Device description	16
3.4	Accessories (optional).....	17
3.5	Functional description	17
3.6	Description of the treatment programmes	18
4	Operation	19
4.1	Safety information for the user	19
4.2	Description of the skin types	20

4.3	Application times	21
4.4	Operating overview	22
4.4.1	Navigation	23
4.5	Functions.....	23
4.6	Start.....	24
4.7	Connect MP3 player.....	25
5	Cleaning and maintenance	26
5.1	Safety information for cleaning and maintenance	26
5.2	Faults.....	27
5.3	Cleaning	27
5.3.1	Disinfection.....	27
5.3.2	Cleaning the surfaces.....	28
5.3.3	Lamps.....	28
5.4	Cleaning plan	29
5.5	Maintenance.....	30
5.5.1	Acrylic glass panel.....	30
5.6	Maintenance plan	31
5.7	Disassembly for maintenance and cleaning work	32
5.7.1	Service kit.....	32
5.7.2	Removing the acrylic glass panel	33
5.7.3	Removing the canopy panel	34
5.7.4	Removing the front panel	36
5.8	Instructions for replacing the lamps.....	37
5.8.1	Collagen lamps.....	37
5.8.2	Starter.....	38
5.9	Cleaning or replacing the lamps in the base	39
5.10	Cleaning or replacing lamps in the canopy.....	40
5.11	Cleaning the base filter.....	41
6	Technical data	43
6.1	Power rating, electrical specifications and noise level.....	43
6.2	Dimensions.....	44
6.3	Lamp sets.....	45
6.3.1	Beauty FLAIR Collagen	45
6.4	Spare parts and accessories	46
7	Appendix	47
7.1	JK timer controls.....	47
7.2	Signs and stickers on the device	48
8	Index.....	51

1 General

1.1 Definitions

Operator

People who provide collagen devices to users on a commercial basis. The operator is responsible for the orderly operation of the device as well as the compliance with the maintenance intervals.

User

A person who uses the collagen device at a commercial location.

Personnel

People who are responsible for the operation, cleaning as well as general maintenance work and who instruct the users in the operation of the devices.

Electrician

A person with suitable professional training, knowledge and experience which enables them to identify and avoid the dangers which can arise from working with electricity.

Authorised, trained and qualified personnel

Qualified staff from an external company who have been trained and authorised by the manufacturer for assembly and maintenance work with specific collagen devices.

1.2 Symbol explanation

The following types of safety notices are employed in these operating instructions:

DANGER!

Type and source of hazard

This safety notice indicates the existence of a direct danger to body and life.

DANGER!

Type and source of hazard

This safety notice warns of dangers to body and life which are caused by electricity.

WARNING!

Type and source of hazard

This safety notice warns of machine, material or environmental damage.

NOTE:

This symbol does not identify any safety notices, but provides additional information to better understand the processes.

2 Information for the operator

2.1 Intended use

This device serves to provide cosmetic stimulation of the skin by means of red light (photo-bio-stimulation) Further information for this can be found on pages 19 and 20.

This device is intended only for commercial use and not for private use.

Some people are not permitted to use the device, see page 19.

The device may not be used by people with reduced physical, sensory or mental capabilities or who have insufficient experience and/or knowledge without supervision. A person responsible for their safety must ensure, either by supervision or instruction, that the device is used correctly and safely. If there is any doubt, use of the device by such categories of persons is to be forbidden!

This device is only to be operated with the specified lamps or their equivalent. The usage times given in these operating instructions are valid only for the specified lamp set.

Any other use shall be considered improper. The manufacturer cannot be held liable for damage or injuries resulting from this. The operator bears the sole risk for this.

2.2 Foreseeable misuse

The following uses of the device are explicitly prohibited:

- Multiple people may not use the device at the same time. The acrylic glass panel may otherwise break and cause injuries.
- People with body weight of more than 135 kg may not use the device. The acrylic glass panel may otherwise break and cause injuries.
- Extending the specified application times cannot place excessive demands on the skin due to the emitted light. Nevertheless we recommend a 48 hour break after a 15-20 minute treatment.
- No other persons are permitted in the vicinity of the cabin when using the device, this especially applies to children.

2.3 Operator's obligations

As the operator you are responsible for providing clear operating, cleaning and maintenance instructions and ensuring the intended use and the proper operation of the device by means of training and instruction for the personnel. Your operating instructions must enable the safe use and safe operation of the device and take into account the characteristics and expertise of your company and also the national work safety and environmental protection regulations. Within the EU the EU directive 89/391/EEC (in Germany the Ordinance on Industrial Safety and Health (BetrSichV)).

2.4 Personnel training

As the operator you are obliged to train your personnel and instruct them in the established legal and accident prevention regulations. Ensure that your personnel have understood and observe these operating instructions. This applies in particular to the

- information in 'Intended use' and 'Foreseeable misuse'
- safety instructions in the 'Operation' chapter
- operating instructions
- cleaning and maintenance instructions
- Please make the instructions and additional information from the manufacturer available to the personnel at their workplace.
- Regularly check your personnel's safety consciousness and awareness of the risks, taking into account your operating instructions.
- Please ensure that the cleaning intervals are complied with, in particular cleaning after every application.
- Please ensure that your staff provide your customers with appropriate advice, in particular regarding application times, possible risks and the operation of the device.

2.5 Scope of delivery

Information regarding the scope of delivery is contained in the 'Description' chapter on page 15.

2.6 Transport, assembly and setup

The device is delivered by a specialist company and assembled and set up by the manufacturer's own service personnel or by an authorised service company.

As an operator you are responsible for complying with the on-site ventilation and electrical requirements at the setup location, see page 9.

2.6.1 Assembly

DANGER!

Personal danger due to electric shock or burns!

- The assembly and electrical connection must comply with the national regulations.
- The on-site electrical installation must be fitted with an easily accessible all-pole isolating device (master switch), complying with overvoltage category III. This means that each pole must have a contact opening width complying with the conditions of overvoltage category III for full isolation.
- If the connection is established via a plug-in connection then the connection system according to EN 60309 must be utilised.
- Assembly, installation, upgrading or repair of the device may only be carried out by specially trained and instructed personnel.

Information regarding the disassembly and assembly of the device is available in the separate assembly instructions. These assembly instructions are intended solely for use by the manufacturer's service personnel or by qualified specialists authorised by the manufacturer.

2.6.2 Ventilation


WARNING!

Risk of overheating due to air flow problems.

Operating faults possible.

- Do not change, obstruct or block the air inflow and outflow to the device or make any unauthorized changes to the device. The manufacturer shall bear no liability for damage or injuries resulting from such actions.
- Observe the ventilation system data.
- Do not set up and operate the device on the transport pallet!

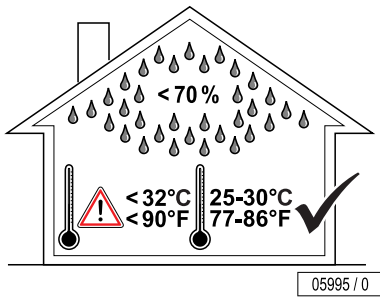
Data regarding the air supply and ventilation planning is available in the Ergoline planning manual (www.ergoline.de, available in English and German). You can also obtain this data from your dealership or the customer service department (see page 2).

2.6.3 Setup location

The device may not be operated outside.

A conversion is required in order to operate the device without problems at altitudes exceeding 2000 meters above sea level. Please contact the customer service department (see page 2) without fail before putting the device into operation.

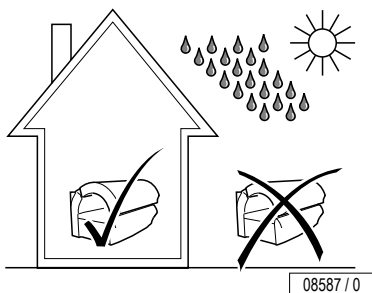2.7 Commissioning

The initial commissioning is performed by the manufacturer's own service personnel or by an authorised service company. The device is handed over ready for use.

If the device has not been used for a longer period of time then it must first be inspected by our customer service department or by another authorised company before being put back into operation again.

The device is operated with an external time-control unit.

A time control unit must be attached so that the device automatically switches off at the latest after <110 % of the selected running time even if the controller fails. The device running time must be provided with double protection by use of a timer control that complies with the standard EN 60335-2-27).

Dependent on the local utilities company, the devices can transmit interference to the power system of the building that impair the audio-frequency remote control system used by the utility company. This can, for example, interfere with the function of night-storage heaters.

If faults occur due to the operation of the devices, then the operator is responsible for the installation of an audio frequency blocking device in the building electrical installation. Please contact your electrician. The electrician is familiar with the technical connection conditions of your local utility company so that the audio-frequency blocking device can be matched to the power system of your utility company.


2.8 Operation and maintenance

Error-free operation, maintenance and service are required in order to ensure that the health and safety of the users is not endangered during operation and to ensure that the device is in a fault-free operating condition.

Please observe the information and recommendations provided in these operating instructions. Please ensure that the intervals for the inspection, maintenance and service are complied with. Please ensure that the lamps are replaced after the specified usage time.

Unauthorized modifications and changes to the device and its control unit are prohibited. Failure to comply with this voids the operating licence!

2.9 Decommissioning

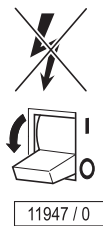The device must be disconnected from the power supply in order to temporarily or permanently decommission it.

The legal requirements regarding disposal must be observed when permanently decommissioning the device.

2.10 Storage

Store the devices in a dry, frost-free location with a stable temperature. Once the device has cooled off, pack it in plastic wrapping to protect it against scratches.

2.11 Disposal

2.11.1 Declaration – The JK Corporate Group

The JK Corporate Group is subject to the strict regulations of EC Directive 761/2001 and the standard EN ISO 14001:2004, and undergoes regular internal and external environment audits performed by trained auditors.

2.11.2 Environmental regulations – disposal of lamps and batteries

Collagen lamps contain fluorescent materials and other waste products containing mercury. Batteries contain heavy metal compounds.

Within the European Union the national implementation of the Waste Framework Directive 75/442/EEG (in Germany the Closed Substance Cycle Waste Management Act (KrWAbgG)) applies. According to the national waste disposal laws and in accordance with the municipal waste regulations, proof must be provided of the proper disposal of lamps and batteries.

Your local sales agency will be happy to assist you with the disposal of lamps and batteries.

- Report the number of lamps and batteries to your local agency by telephone or in writing.

The agency will then provide you with the name of a free delivery point for your lamps or will take care, in conjunction with a disposal company, of the collection and proper disposal of the lamps. This service is available for a fee. Outside Germany, the respective national legislation applies. Please contact your local sales agency.

2.11.3 Packaging

All packaging consists of 100 % recyclable materials. Packaging brought into circulation by the JK Corporate Group that is no longer required can be returned to the JK Corporate Group. Your partner agency or dealer will be happy to advise.


2.11.4 Disposal of old devices

The device has been manufactured using recyclable materials. The JK Corporate Group will provide you with information on the content or potential hazards of the materials used.

In accordance with the 2002/96/EC directive, implemented in Germany in the Electrical and Electronic Equipment Act (ElektroG), the manufacturer is obliged to take back and dispose of certain electrical and electronic components.

The components and devices are labelled with the following symbol:

The device will, if requested, be properly disposed of by the JK Corporate Group. This service is available free of charge. Your partner agency or dealer will be happy to advise.

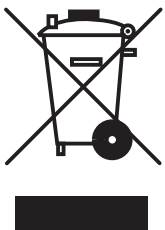2.12 Directives

This device was built according to the following directives:

- EC directive 'Electromagnetic Compatibility' 204/108/EC (according to the version valid at the time).
- Low voltage directive 2006/95/EC (according to the version valid at the time).

2.13 Export

We emphasize that these devices are intended for the European market only and must not be exported to or operated in the USA or Canada! The manufacturer does not accept any liability in case of non-observance! It is emphatically emphasized that non-compliance with this instruction may result in high liability risks for the exporter and/or the operator.

2.14 MP3 music (optional)

If an MP3 player is used privately during the application, then this does not constitute a public performance in the sense of copyright law, so that no reporting or fee obligation exists on the part of the studio operator towards the organisation responsible for issuing the commercial performance rights.

For the public use of MP3 music files, the same legislation applies as for all other music sources: As a studio operator/studio owner, you may only play original Audio-CDs, MCs, audio-DVDs etc., on your premises and/or play music via the 'MP3-music modules' which are provided in the professional devices in the Ergoline range, if you have acquired the necessary stage rights. You can acquire these from the organisation responsible for issuing commercial performance rights in your country (in Germany: GEMA/GVL).

Copy-protected audio-CDs, MCs, audio-DVDs etc., as well as titles contained on these media, may not be converted to MP3 format and/or saved on HDD, audio-CDs, MCs, audio-DVDs etc., if software is used to do this, which disables or circumvents the copy-protection contained on the media.

You must be ready at any time to prove your granting of stage rights by presentation of the appropriate stage rights certificate when requested by a representative or another controlling body of the responsible organisation.

As long as you observe all the listed instructions and have acquired the necessary rights for use of the 'MP3-music-module', you may convert your original audio-CDs, MCs, audio-DVDs etc. into the MP3 format required by the player. Only one copy (duplicate or MP3 format conversion) of lawfully acquired audio-CDs, MCs, audio-DVDs etc. may be made. The source media (original media) must be retained and may not be used simultaneously.

Further details can be found in the information brochure 'Audio Guide' available in German and English on the Internet at www.jk-global-service.de.

2.15 Instruction signs on the device

See Appendix page 48.

2.16 Technical modifications

The device has been constructed in accordance with the current state-of-the-art and the applicable safety regulations. The illustrations and specifications on these instructions are subject to technical modifications which are required in order to improve the device.

Unauthorised modifications and changes to the device and its control unit are prohibited. Failure to comply voids the EC conformity declaration and the GS seal!

Please use only original spare parts. Failure to comply voids the EC conformity declaration and the GS seal!

2.17 Manufacturer's guarantee

On behalf of customers who have purchased an Ergoline collagen device for private or commercial utilisation from one of Ergoline's sales partners, Ergoline accepts liability for any defects that arise in accordance with the following regulations; consumable parts such as collagen lamps and starters as well as acrylic glass panels are not covered by the warranty.

The guarantee stipulates that Ergoline will rectify any defects within a time period predetermined by Ergoline by repair or exchange of the defective part. Claims made in accordance with this guarantee can be lodged with Ergoline or the Ergoline sales partner from whom the product was purchased on presentation of the guarantee card or invoice; this guarantee is valid for a period of 24 months from the date of purchase.

This does not affect the rights of the customer under the existing guarantee with regard to any claims that the customer might have in accordance with the contract of purchase.

2.18 Warranty and liability exclusions

Warranty and liability claims for personal injury and property damage are excluded if they are the result of one or more of the following causes:

- improper use of the device;
- improper assembly, commissioning, operation and maintenance of the device;
- operating the device with faulty safety equipment or improperly attached or non-functional safety and protective equipment incl. warning signs;
- not complying with the information in the operating instructions regarding transport, storage, assembly, commissioning, operation and maintenance;
- the use of untrained personnel;
- unauthorised changes to the device or the control unit;
- inadequate monitoring of parts subject to wear;
- repairs carried out improperly;
- the use of non-original spare parts;
- emergencies and disasters resulting from the influence of foreign materials or force majeure.

3 Description

3.1 Scope of delivery

- Collagen device
- Technical documentation (folder with instructions, brochures for error codes and default settings and additional documentation)
- Service kit: with allen key and suction cup for canopy panel and acrylic glass panel
- Collagen lamps
- Connecting lead

3.1.1 Optional

- Stereo Sound System or Stereo Sound System PLUS
- MP3 cable


3.2 Equipment

Device equipment varies depending on the model.

In this description, all standard and optional components/functions are taken into account which must be observed during operation and/or when performing maintenance.


The control panel also differs according to the varying equipment: Only those buttons are visible that can be operated, i.e. the corresponding function must be included in the device's equipment.

3.3 Device description

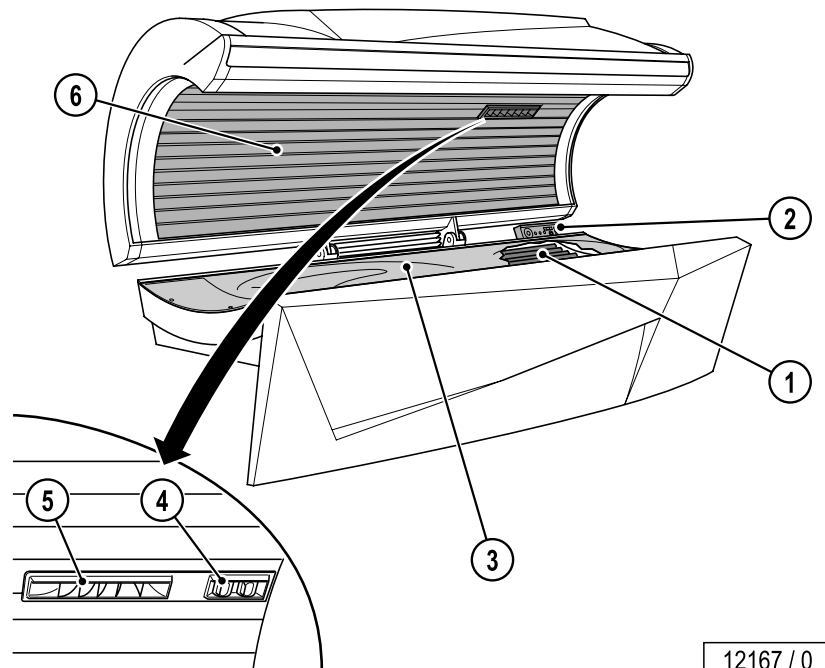1. Collagen lamps, base
2. Operating panel
3. Acrylic glass panel base
4. Adjustable air nozzles face cooling
5. Air nozzles body cooling
6. Collagen lamps, canopy

3.4 Accessories (optional)

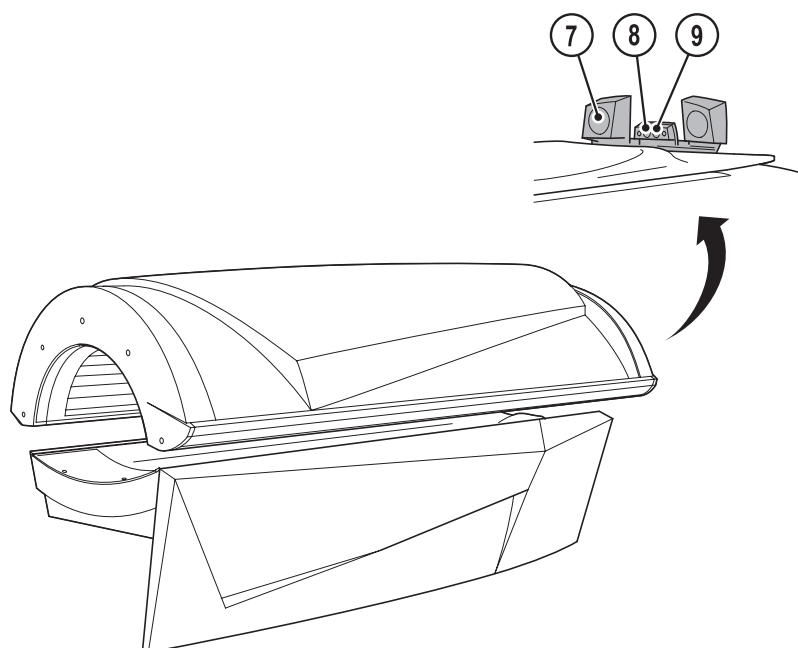12184 / 0

- 7. Audio system (Stereo-Sound system PLUS or Stereo-Sound system)
- 8. Headphone jack (only with Stereo-Sound system PLUS)
- 9. MP3 connection (only with Stereo-Sound system PLUS)

3.5 Functional description

The application with red light (photo-bio-stimulation) promotes the formation of collagen and elastin.

The short-term effect shows itself as a firmer and closer skin feeling and a more vital facial colouration. This is the result of oxygen enrichment and an increased level of ATP (adenosine triphosphate, the molecule for energy transfer in the metabolic system). During the application process the serotonin level also increases which results in a feeling of wellness and relaxation.

Over the middle and long-term the skin complexion improves as a result of greater moisture and oxygen levels. Natural signs of ageing are mitigated, the skin's collagen and elastin network strengthened and revitalised.

A complete treatment consists of at least 8 - 12 applications over a period of approximately 4 weeks. A higher number of applications over a longer period of time is recommended for more mature users.

Following the treatment 1 - 2 applications per month are sufficient to maintain the effect. The effect is supported by the use of anti-ageing skin care products.

3.6 Description of the treatment programmes

AGEsmart programme RENEW	AGEsmart programme VITALIZE	AGEsmart programme PROTECT
Suitable for more mature users in order to alleviate wrinkles and actively combat skin aging	Suitable for mature users who wish to refresh their skin and counteract the initial signs of aging	Ideal for younger users of 18 years and over for refining their complexion and protecting against environmental influences
1. Stimulates the collagen and elastin production and vitalises the skin	►	►
2. Improves the absorption and effect of cosmetic products	►	►
3. Refines the pores and improves the complexion	►	►
4. Provides the skin with a more lustrous and healthy appearance	►	►
5. Increases the resistance to environmental and stress-related influences	►	►
6. Detoxifies the skin and encourages cell regeneration	►	–
7. Increases the skin's oxygen and moisture levels	►	–
8. Reduces pigmentation and liver spots	►	–
9. Reduces swelling and dark rings around the eyes.	►	–
10. Substantial reduction of the wrinkles in facial and neckline areas	Counteracts fine facial lines and wrinkles	–
11. Fills the connecting tissues with collagen and elastin – for visibly firmer skin on the entire body	–	–
12. Effective and long-term reduction of skin aging	–	–

NOTE:

The application times can be found on page 21.

4 Operation

4.1 Safety information for the user

DANGER!

Danger of hearing damage!

Using earphones or headphones at high volume can result in permanent hearing damage.

- Do not set the volume too high.

NOTE:

Extending the specified application times cannot place excessive demands on the skin due to the emitted light. However we recommend a 48 hour break after a 15-20 minute treatment.

The following people are not allowed to use the equipment:

- Pregnant women
 - People younger than 18
 - People with epilepsy
 - People with illnesses which cause light sensitivity
 - People who take medication or herbal remedies which cause light sensitivity (e.g. antibiotics, antirheumatic agents, antifungal agents or St. Johns wort)
- **Seek medical advice immediately if in doubt!**

A light treatment is not recommended in the following circumstances:

- If you have utilised an oxygen facial pack in the last 24 hours,
- If your facial skin has been electronically tightened in the last 24 hours,
- If your skin has been scoured,
- If you have received a facial injection in the last 24 hours.

Behaviour before, during and after a light application:

- No other persons are permitted in the vicinity of the cabin when using the device, this especially applies to children.
- Remove all jewellery (including piercings) prior to the application.
- Remove cosmetics before using the device and do not use sun screen products.
- A light application is not damaging to the eyes. However, we nevertheless recommend wearing eye protection during the application in order to protect the eyes against the bright light.
- Always comply with the recommended application times.

4.2 Description of the skin types

NOTE:




Extending the specified application times cannot place excessive demands on the skin due to the emitted light. However we recommend a 48 hour break after a 15-20 minute treatment.

- Comply with the application times.
- Observe the safety instructions.

NOTE:

The application times are only valid for the lamp set specified on the sticker.

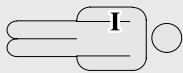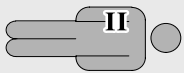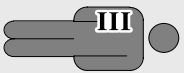
Skin type I	Skin type II	Skin type III	Skin type IV
			
Sensitive skin	Light skin	Normal skin	Dark skin
Always/frequently sunburn	Frequently sunburn	Rarely sunburn	Rarely sunburn
Tolerates only a small amount of natural sun.	Tolerates approx. 10-20 min. of natural sun.	Tolerates approx. 20-20 min. of natural sun.	Tolerates approx. 40 min. of natural sun.

4.3 Application times


AGEsmart programme –Treatment recommendations

Applications per week	PROTECT		VITALIZE		RENEW	
	Light	Intensive	LightLight	Intensive	Light	Intensive
	15-20 min	15-20 min	15-20 min	15-20 min	15-20 min	15-20 min
	[Number]	[Number]	[Number]	[Number]	[Number]	[Number]
1. Week	2	2	2	2	2	2
2. Week	2	2	2	2	2	2
3. Week	2	2	2	2	2	2
4. Week	2	2	2	2	2	2
5. Week		1	2	2	2	2
6. Week		1	2	2	2	2
7. Week				1	2	2
8. Week				1	2	2
9. Week					1	2
10. Week					1	2
11. Week						1
12. Week						1
After treatment:	1-2 x per month		1-2 x per month		1-2 x per month	

4.4 Operating overview

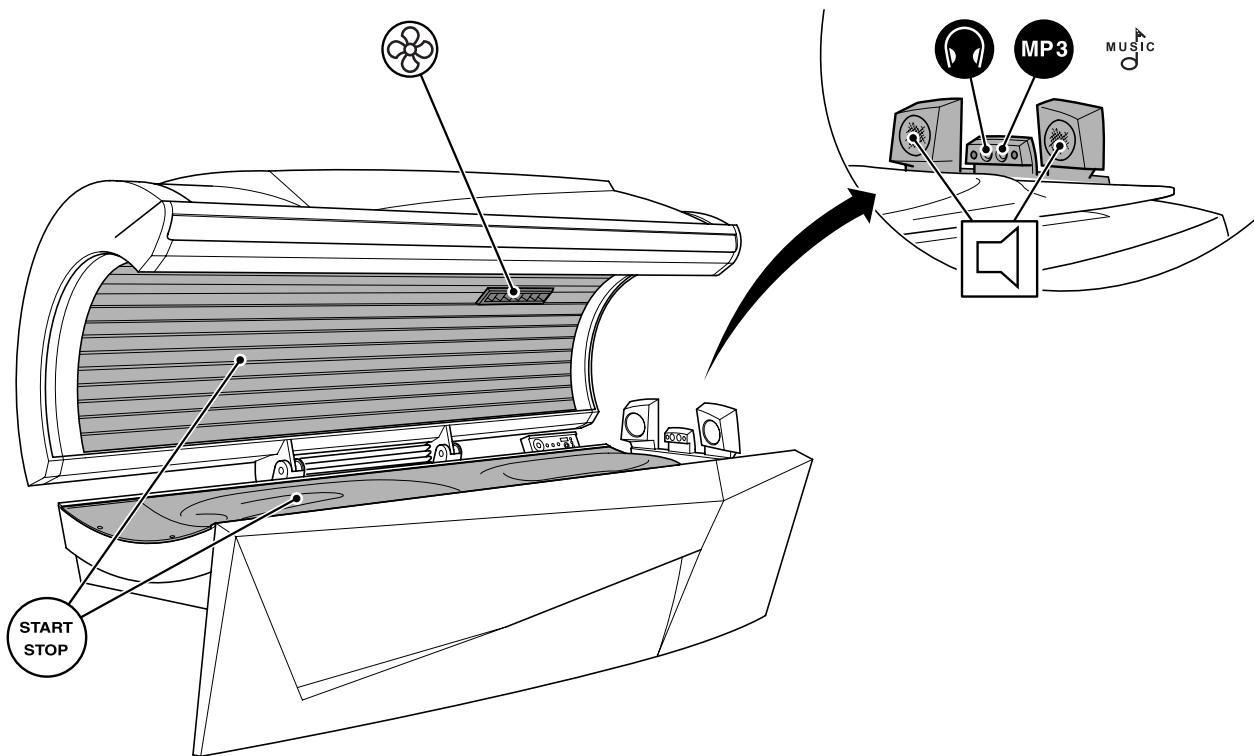
12169 / 0

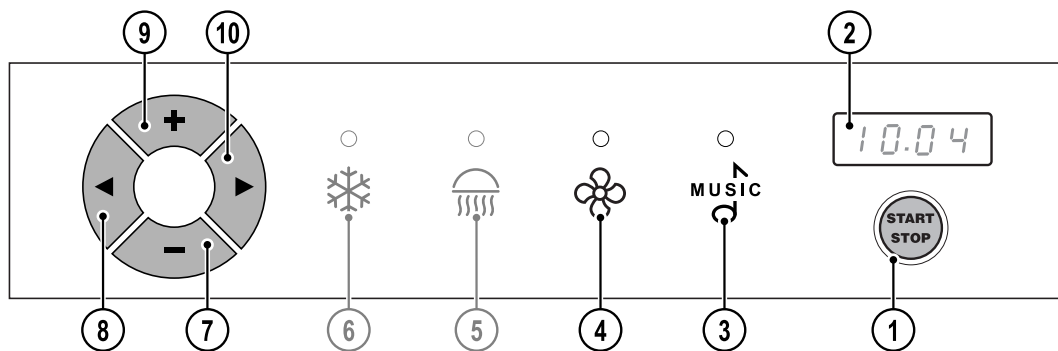12238 / 0

1. START/STOP button
2. Display – Functions ...
3. Music pictogram (only available if Stereo Sound or Stereo Sound PLUS is present)
4. Body cooling/facial cooling pictogram
5. Facial tanner pictogram (not available)
6. Air conditioning pictogram (not available)
7. Minus button
8. Left navigation arrow button
9. Plus button
10. Right navigation arrow button

4.4.1 Navigation

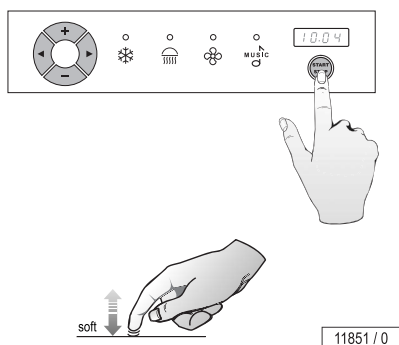The device is operated with the navigation ring.

Switch on lamps:

- Tip the button once.


Switch off lamps/interrupt application:

- Hold the button.

Select device function (go to pictogram)

MUSIC e.g. MUSIC

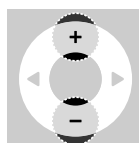Select sub-function (if available)


Change setting

The diodes of the available functions and the navigation ring are illuminated in default mode after starting up the device. After selecting a function, the respective diodes and the navigation ring flash for a few seconds. Settings can be changed during this time.

4.5 Functions

Function	Button sequence	Description	Information
1. Lamps		START/STOP Switch the lamps on / off during the application	If the lamps are switched off during an application, the time continues to run.
3. Music		Internal: Control volume, Levels 0-14 External: MP3 player	Only available with optional accessories Stereo-Sound and Stereo-Sound PLUS. See page 25
4. Body cooling/ Facial cooling	 	Control cooling, Levels 1-5 Switch off cooling <input type="text" value="OFF"/> Switch on cooling <input type="text" value="ON"/>	The air jets for the facial cooling can be adjusted by hand.

4.6 Start

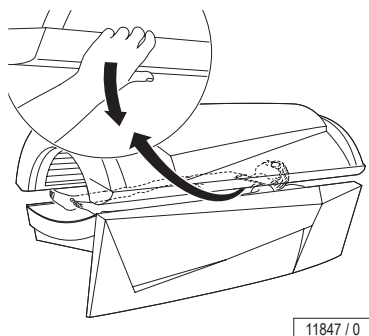- Close the canopy.
- Tip the START/STOP button.

Upon starting the following functions are switched on:

- Lamps
- Body cooling
- Stereo-Sound system PLUS (optional)

4.7 Connect MP3 player

DANGER!

Danger of hearing damage!

Using earphones or headphones at high volume can result in permanent hearing damage.

- Do not set the volume too high.

Users may connect their private MP3 players to the device if the 'Stereo-Sound system PLUS' accessory is fitted.

Mat.-Nr. 1001855-... provided for use by the studio operator
(1.5 m / 2 x 3.5 mm stereo jack plug).




MP3 player connection: see the documentation from the manufacturer of the MP3 player.

Connect:

- Plug the MP3 cable (stereo jack plug) into the Stereo-Sound system PLUS.
- Plug the MP3 cable into the MP3 player.
- Select the music title and start.

If the speakers do not play any music:

- Increase the volume on the MP3 player until the music signal is automatically switched to the Stereo-Sound system PLUS.
- Set the volume on the device to the desired level.

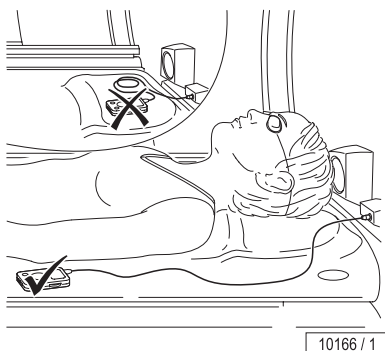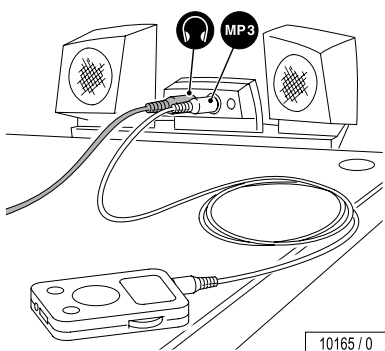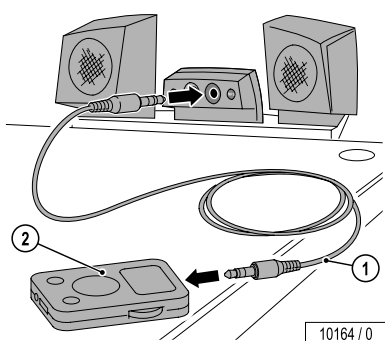WARNING!

Heat damage to the MP3 player is possible

- Do not position it directly in the radiation area.

Cable may break-off, cable jacks may be damaged!

- Do not pull on the cable when unplugging the MP3 player.

5 Cleaning and maintenance

5.1 Safety information for cleaning and maintenance


DANGER!

Electricity throughout entire device

Personnel danger due to electric shock or electric burns.

- Before beginning work disconnect all of the power connections from mains voltage.
- Secure all of the disconnected connections against accidental reactivation.
- Check that the device is disconnected from mains voltage.
- All of the work and electrical connections must comply with the national regulations of the country in which the device is being operated (e.g. the VDE regulations in Germany) and must be performed by appropriately trained specialists.

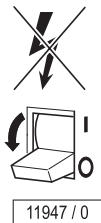
If work is to be performed on the devices, they must be deenergized. That means that all live cables must be switched off.

Switching off the device alone is insufficient, as certain points may remain energised. Therefore, switch off fuses before performing any work and – if possible – remove them.

Accidental reactivation can result in serious accidents. Immediately after deenergising, secure all switches or fuses used for deenergising against being switched on again.

- Lock the fuse box with a padlock.

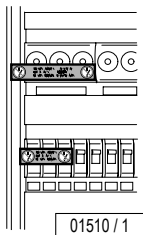In case of circuit breakers which cannot be unscrewed, a strip of adhesive tape may be applied over the actuation lever with the wording 'Do not switch, danger!'.

A prohibition sign with the text:

'Work in progress!'

'Place:'

'Sign may only be removed by:'

must always be securely mounted immediately.

5.2 Faults

Error codes are output to the display to simplify troubleshooting:

- If a fault occurs, the error code flashes on the display.
- If several errors occur, the error messages are displayed alternately.
- Error elimination is acknowledged by pressing the **START/STOP** button.
- If the error cannot be cleared, inform Customer Service - see page 2.

More detailed information regarding the error description and solutions is available in 'Fehlercodes/Error codes', (Order-no. 1008362-..).

5.3 Cleaning

DANGER!

Infections can be spread by skin contact!

All objects/device parts, which can come into contact with the user during the application, must be disinfected after every application:

- Acrylic glass panel
- Handles and control panel
- Protective goggles
- MP3 cable

5.3.1 Disinfection

We recommend the following cleaning agents for rapid and thorough disinfection:

Antifect® fast disinfectant cleaner

Antifect® concentrate, 250 ml	Order-no.:	8001813-..
Spray bottle, 1 ltr. (empty)	Order-no.:	8001513-..
Spray nozzle	Order-no.:	8001613-..
Mixing container, 5 ltr. (empty)	Order-no.:	8001713-..

NOTE:

For thorough disinfection, the application time must be adhered to. Follow the manufacturer's instructions for use.

5.3.2 Cleaning the surfaces

Acrylic glass panel

WARNING!

Do not rub with a dry cloth – danger of scratching!

Failure to comply voids any warranty claims.

- For fast, hygienic proper cleaning of acrylic glass surfaces, use the specially developed fast disinfectant cleaner Antifect®.
- Other cleaning agents, especially concentrated disinfectants or solvents (e.g. Lysoform, ethyl alcohol or other liquids containing alcohol) must never be used for cleaning.

Plastic surfaces

For cleaning the remaining painted and unpainted plastic surfaces, it is best to simply use warm water and a leather cloth. Never use aggressive cleaning agents containing alcohol or essential oils. These may lead, in the long run, to damages not covered by the guarantee.

During cleaning a slight discolouration of the leather cloth by the rubber seals must be expected for production-related reasons.

Prevent damage to the acrylic glass and plastic surfaces.

- Remove rings, watches, bracelets etc. before starting cleaning work.


02756 / 0

5.3.3 Lamps

Lamps can be cleaned with a damp cloth and without additional cleaning agents where necessary.

5.4 Cleaning plan

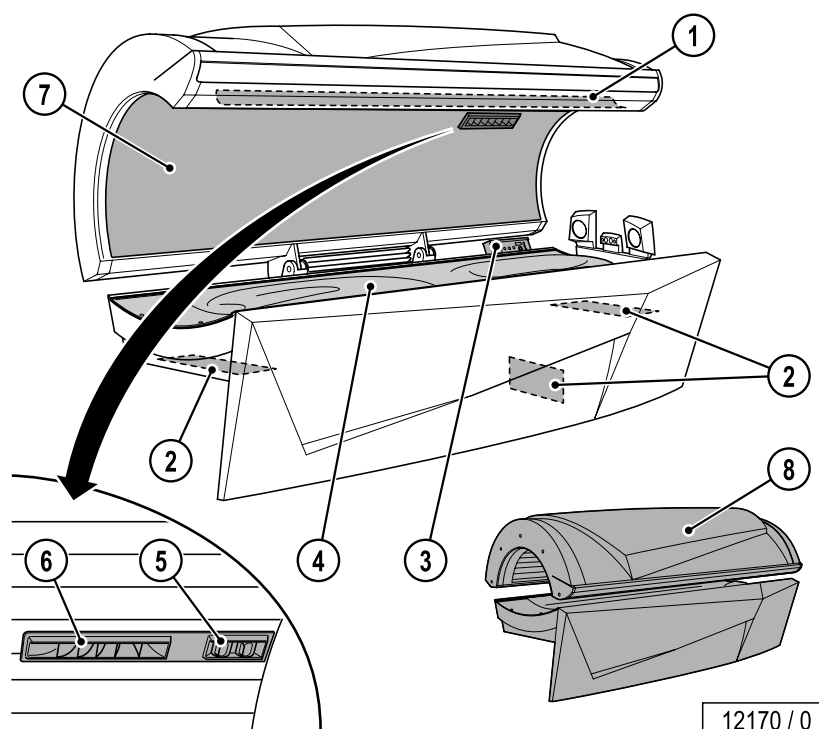Component	Interval	Activity	See page
	After each application When necessary 50 h		
1 Handle	X	Cleaning and disinfection	27
3 Control panel	X	Cleaning and disinfection	27
4 Acrylic glass panel	X	Cleaning and disinfection	27
5 Adjustable air nozzles	X	Cleaning and disinfection	27
Protective goggles	X	Cleaning and disinfection	27
MP3 cable	X	Cleaning and disinfection	27
6 Air jets	X	Clean	27
7–8 Surfaces	X	Clean	27
2 Filter	X	Remove and clean filter	41

NOTE:

Cleaning must be confirmed after each application by pressing the START/STOP button (when the warm-down time for the main fan has finished).

5.5 Maintenance

Maintenance and care play a major role in deciding whether the device meets the requirements placed on it in a satisfactory manner. Compliance with the specified maintenance intervals and the careful conducting of the maintenance and care work are therefore absolutely necessary.

The current operating hours of maintenance parts can be called up in service mode (see 'Voreinstellungen/Presettings', Order-No. 1008368-..).

We point out that the device must be subjected to regular inspections by our customer service or an approved specialized company every 12 months (from commissioning) to maintain its proper condition!

5.5.1 Acrylic glass panel

Over the course of time the acrylic glass panel becomes less translucent as a result of usage and must be replaced on regular basis.


5.6 Maintenance plan

WARNING!

Only use genuine **spare parts** of the same type! If other parts are used then the device no longer complies with the conformity declaration and may not be operated any longer!

Any liability for damage or injury proven to be the result of the use of non-genuine spare parts is excluded.

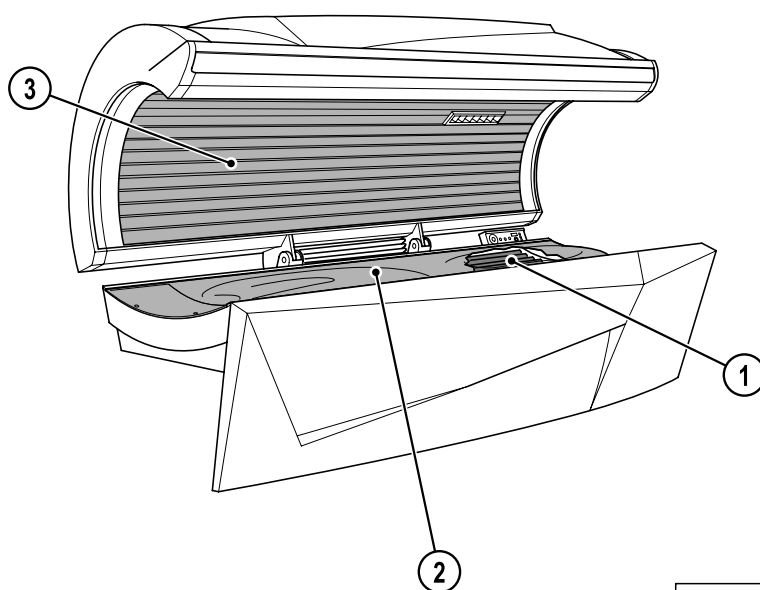12171 / 0



Component	Interval			Activity	See page
	500 h	1000 h	1500 h		
1 Collagen lamp base	X			Replace	39
3 Collagen lamps canopy	X			Replace	40
1 Starter		X		Replace	39
3 Starter		X		Replace	40
2 Acrylic glass panel			X	Replace	32

5.7 Disassembly for maintenance and cleaning work


When necessary the acrylic glass panel, the canopy panel or the front panel can be removed

5.7.1 Service kit


1. Allen key
2. Suction cup for acrylic glass panel and canopy panel
3. Scratch prevention for the entry area (support cushion for the canopy)

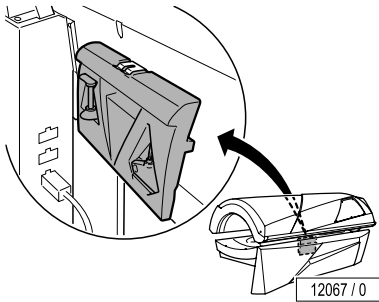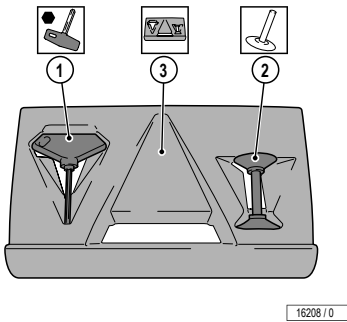
- The service kit is fitted to the rear wall of the device.

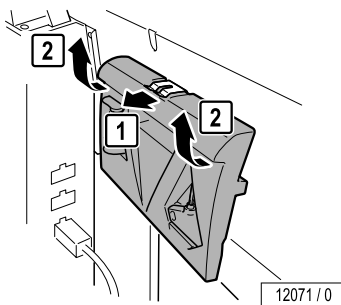- Remove service kit from support.

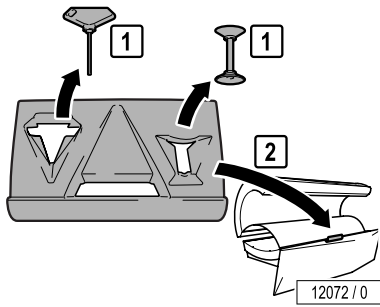
- Remove Allen key and suction cup.

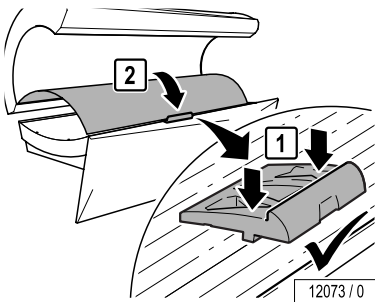- Use cushion as scratch prevention when opening the canopy panel.


5.7.2 Removing the acrylic glass panel

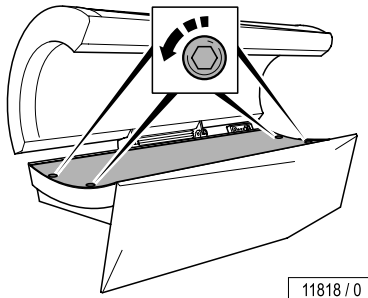
- Remove the screws and remove the washers.

- Push the acrylic glass panel to the rear.
- Use the suction cup to lift the panel out of the guide at the front.

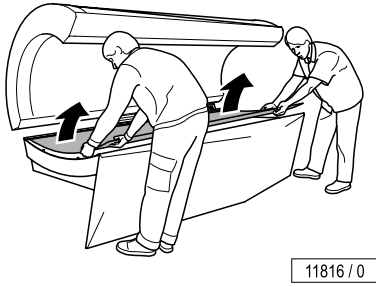- Remove the acrylic glass panel and carefully lay it aside.

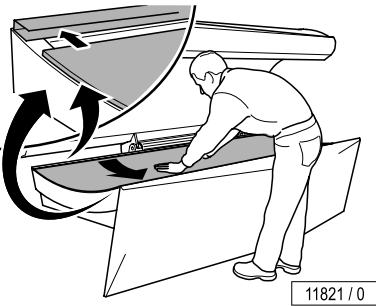
The assembly is the same process in reverse order.


When assembling the sunbed please note:

- Carefully lay the acrylic glass panel in the guide.


5.7.3 Removing the canopy panel

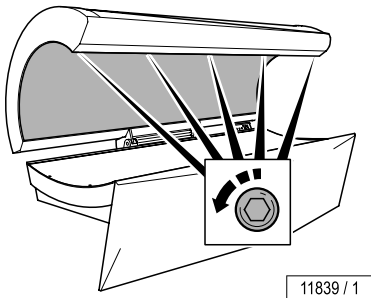
- Remove the screws. First remove the outside screws and then remove the screws in the middle.

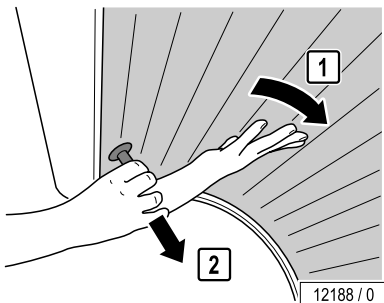
- Push the canopy panel backwards with one hand and use the suction cup to lift it out of the front guide.

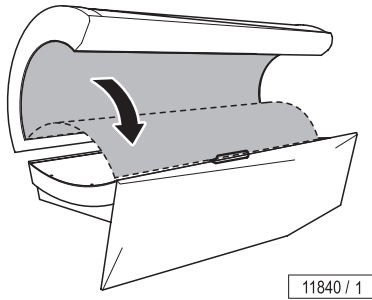- Carefully lay the panel aside on the support cushion (service kit).
- Protect the acrylic glass panel against scratches.

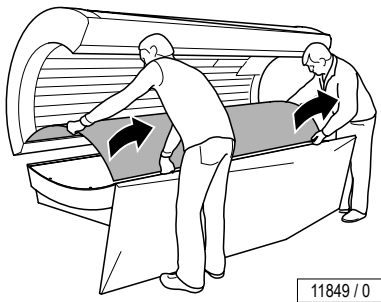- Remove the canopy panel with the aid of a helper.

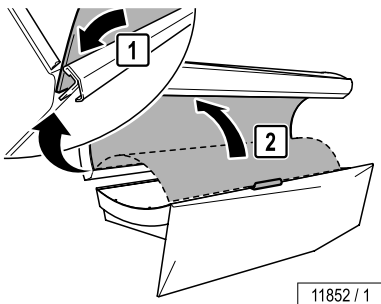
The assembly is the same process in reverse order.


When assembling the sunbed please note:

- Use scratch protection.
- Carefully lay the canopy panel in the guide.


5.7.4 Removing the front panel

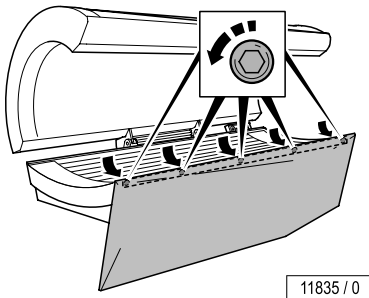- Remove the acrylic glass panel: see page 33.
- If necessary, loosen the screws from the inside.

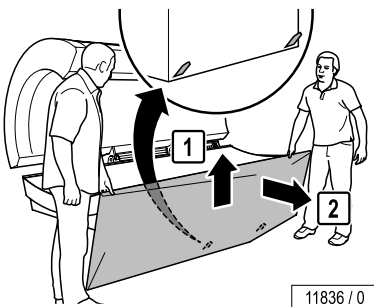- Carefully lift the upper edge of the front panel out of clamps (pull upwards).
- Lift the front panel in order to lift it out of the two lower hooks.
- Then remove the front panel to the front and lay it aside.

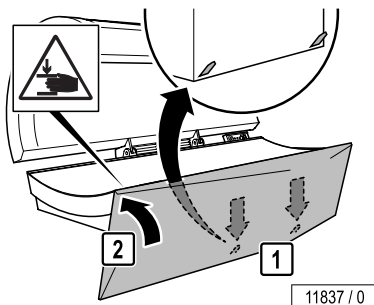The assembly is the same process in reverse order.

DANGER!

Danger of crushing fingers when inserting the front panel!

Fingers can be jammed when pressing in the front panel.



- Do not hold fingers between the front panel and the device when pressing.

5.8 Instructions for replacing the lamps


The effect of the lamps decreases over the course of time. Therefore, the replacement intervals specified in the maintenance plan must be adhered to.

5.8.1 Collagen lamps

- Replace faulty lamps immediately.
- Always replace the entire set of lamps after the specified operating time.
- Turn 90° and carefully remove from the socket.

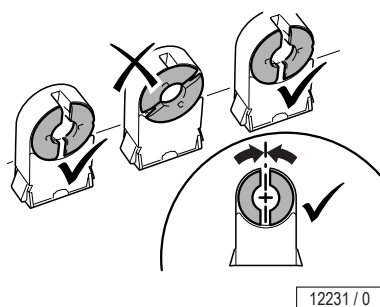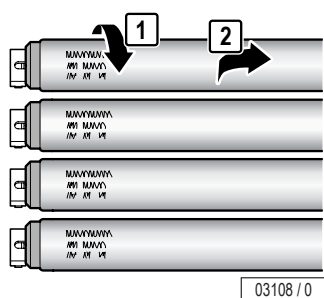- Push the pins vertically into the socket from above.

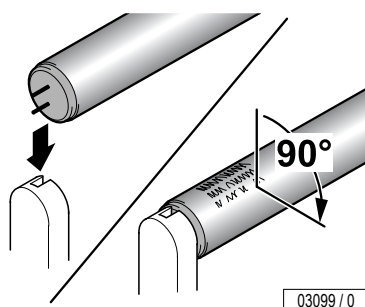- Turn the lamp 90°. The lettering must be oriented towards the outside (i.e. towards the user).


5.8.2 Starter

- Replace faulty starters immediately.

You can recognise a faulty starter by the blackened ends of the UV low-pressure lamps or when the lamp can no longer be started.

Operating a lamp with a faulty starter significantly shortens its service life. The lamp loses power and may begin to flicker under certain circumstances.


The starters are located in either one of the two lamp sockets or are visible when the lamps are removed.

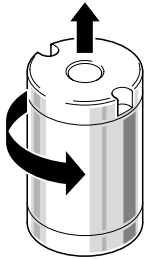05584 / 0


- Turn 90° and carefully remove from the socket.

5.9 Cleaning or replacing the lamps in the base

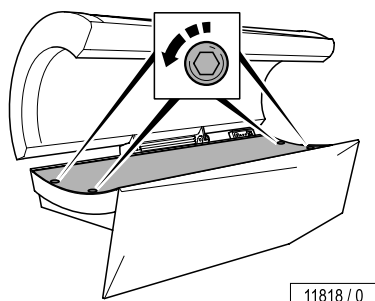11818 / 0

- Remove the acrylic glass panel: see page 33.

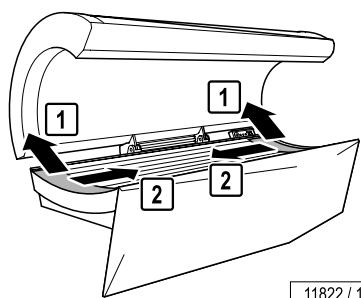
11822 / 1

- Remove the lamp coverings at the head and foot ends.

Replace the lamps: see page 37.

Replace the starter: see page 38.

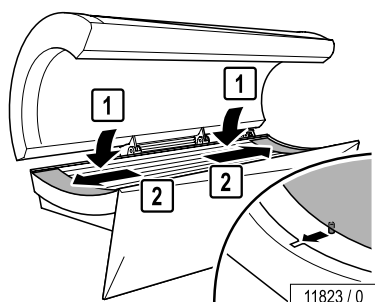
11823 / 0

The assembly is the same process in reverse order.

When assembling the sunbed please note:


- When positioning the lamp covering push the pins into the opening.
- Carefully lay the acrylic glass panel in the guide.

5.10 Cleaning or replacing lamps in the canopy

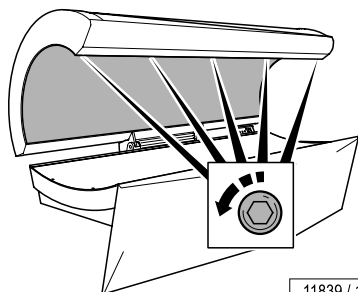
11839 / 1

- Remove the canopy panel: see page 34.

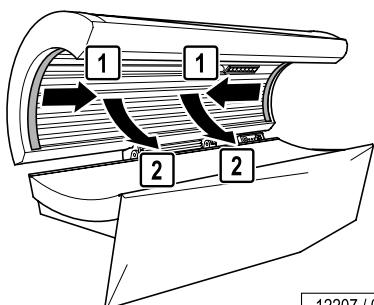12207 / 0

- Remove the lamp cover.
- Replace the lamps: see page 37.
Replace the starter: see page 38.

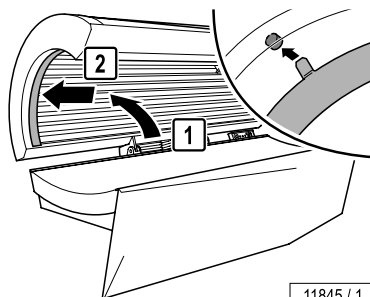
11845 / 1

The assembly is the same process in reverse order.

When assembling the sunbed please note:


- When setting the lamp covers, ensure that the tabs are in contact with the magnetic catch.

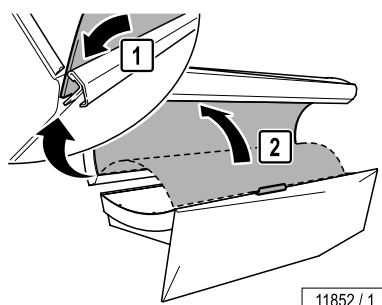11852 / 1


- Carefully lay the canopy panel in the guide.

5.11 Cleaning the base filter

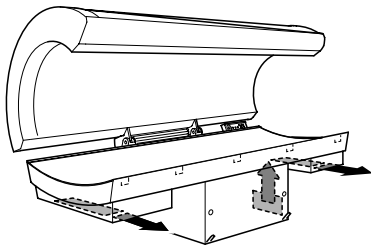
11838 / 0

- Remove the front panel: see page **Fehler! Textmarke nicht definiert..**
- Pull out the two outside filters to the front. Remove the middle filter via the top

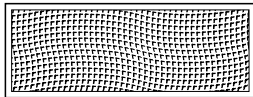
05064 / 0

- Moist or dry clean the filter.
Dry cleaning: Vacuum cleaner (depending on the degree of soiling)
- Moist cleaning: Water and detergent, or with a dishwasher.

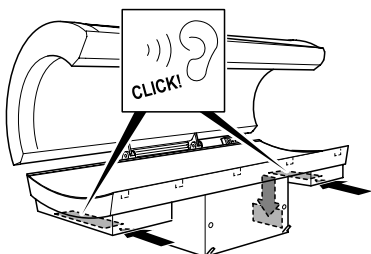12199 / 0

- Replace the cleaned and dried filter. The two outer filters click into place audibly.

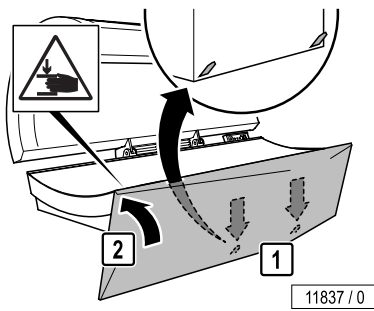- Replace the front panel.

DANGER!

Danger of crushing fingers when inserting the front panel!

Fingers can be jammed when pressing in the front panel.

- Do not hold fingers between the front panel and the device when pressing.

WARNING!

Damage to the device through moisture possible.

The cleaned filters must be dry before they are replaced

6 Technical data

6.1 Power rating, electrical specifications and noise level

The Safety Checked seal (GS seal) certifies that the device fulfils the requirements of the Equipment and Product Safety Act (GPSG). The test mark is located on the type plate.

Beauty FLAIR Collagen

Nominal power consumption without air conditioning:	3300 W
Rated frequency:	50 Hz
Rated voltage:	230-240 V ~1N
Rated fuses:	1 x 16 A (delayed)
Connecting lead:	safety plug (type F) ¹⁾
Noise level at 1 m distance from the device:	62,4 dB(A)
Noise level in device:	<81 dB(A)
Load on the acrylic glass panel:	max. 135 kg
Weight:	310 kg




- 1) Also supplied is a power supply cable with a earth contact type plug (type F).

The following countries will require an adapter to convert a type F to a type G or J plug:


Great Britain: 2 flat contacts and 1 rectangular contact (type G) with 6.3 A micro-fuse

Switzerland: Hexagonal plug with 3 round contacts arranged in a triangular pattern (type J)

Austria: no adapter necessary (type F)

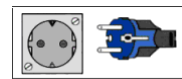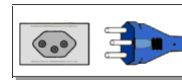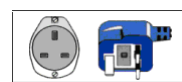6.2 Dimensions

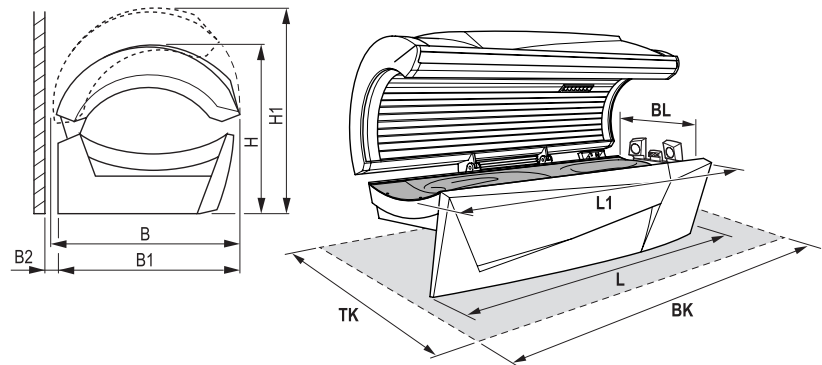12172 / 0

B =	1100 mm
B1 =	710 mm
H =	1100 mm
H1 =	1470 mm
L =	2150 mm (without Sound system)
L =	2250 mm (with Sound system)
L1 =	2150 mm
BK =	2400 mm
TK =	2300 mm


6.3 Lamp sets

WARNING!

The conformity declaration and the GS test mark only apply for the specified lamps. The use of other lamps voids the operating licence.


6.3.1 Beauty FLAIR Collagen

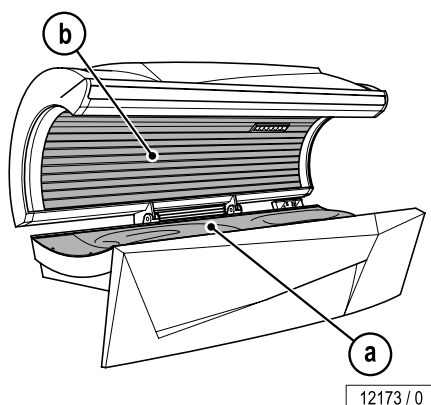Lamp set

	Number	Description	Length	Power	Control	Order-no.
a)	12	Ergoline COLLAGEN BEAUTY LIGHT	1,80 m	100 W	100 W	1504742-..
	12	Starter Q12				1004929-..
b)	16	Ergoline COLLAGEN BEAUTY LIGHT	2,00 m	100 W	100 W	1504741-..
	16	Starter Q12				1004929-..

6.4 Spare parts and accessories

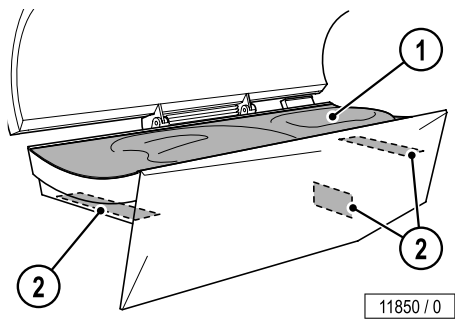Number		Description	Order-no.
1	1	Acrylic glass panel	1503186-..
2	3	Filter	55567-..

7 Appendix

7.1 JK timer controls

The following coin modules/time control systems can be utilised with the Beauty FLAIR series devices.

MCS III plus
Order-No.: 500000456

MCS IV plus
Order-No.: 34010400

MCS VI plus
Order-No.: 34009700


Studiopilot
Order-No.: 34009900


7.2 Signs and stickers on the device

Signs warning of danger areas are attached to the device. The signs shown below are examples. Ensure that the instruction signs are always clearly recognizable and easily legible.

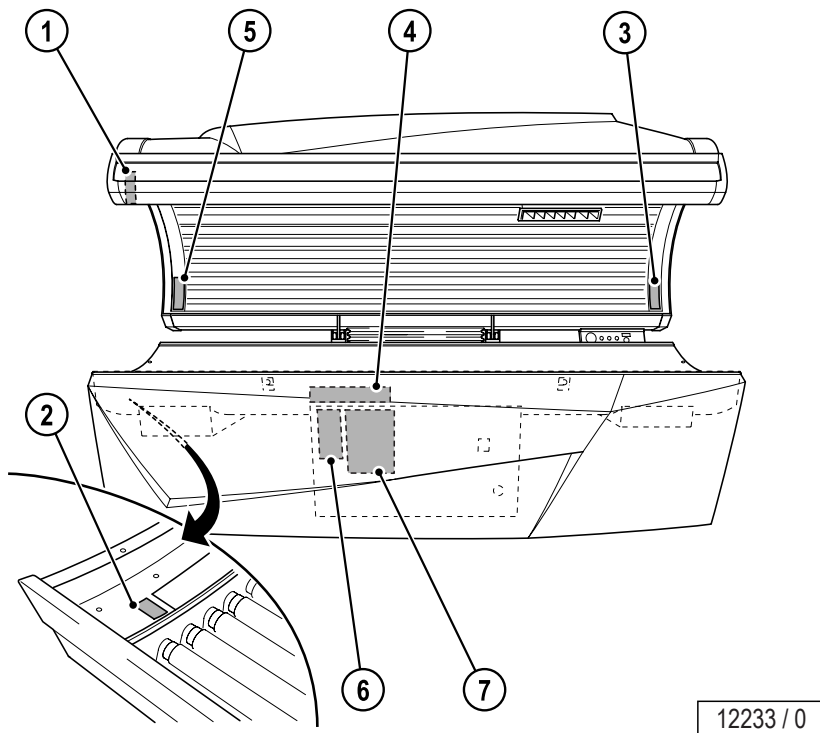1 Name plate

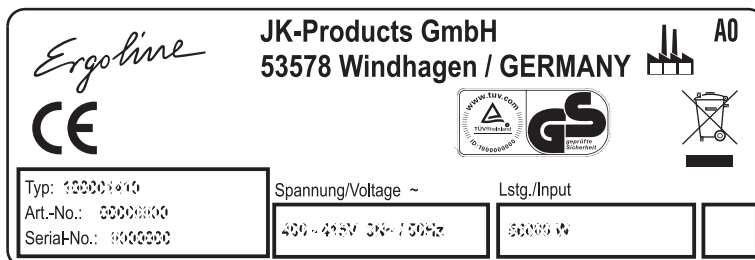The name plate is attached to the interior side of the canopy. It contains important information for identifying the device (e. g. the serial number = Serial-No.).

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12234 / 0

The sticker is located below the lamp covering at the foot end.

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NOTE

Please note that these sunbeds must not be exported to the USA* or Canada or operated there. No liability will be accepted in the event of any failure to observe this notice.

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
84829-02-

11977 / 0

The sticker is located behind the front panel.

5 Maintenance instruction sticker (800701-..)

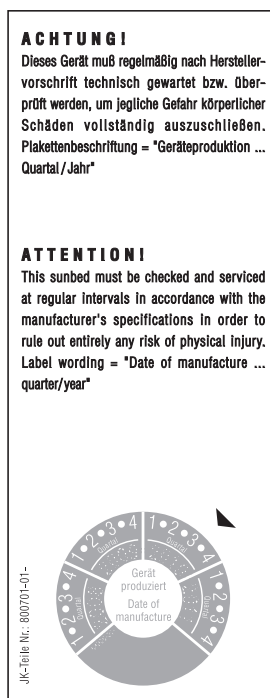12225 / 0

6 Mains power connection sticker 230/400 V devices (1504659-..)

The sticker is visible when the front panel is removed and the control unit cover is opened.

7 EU power supply units sticker

The sticker is visible when the front panel is removed and the control unit cover is opened.

8 Index

A

Accessories.....	17, 46
Acrylic glass (maintenance)	28
Acrylic glass panel	
Permitted load	43
Acrylic glass panel (product information).....	30
Adolescents	7, 19
AGEsmart programmes	18
Application times	21
Assembly	9

C

Cleaning.....	26, 27
Filter	41
Lamps	28
Cleaning intervals	29
Cleaning or replacing the lamps	
Base	39
Canopy	40
Cleaning overview.....	29
Collagen lamps, maintenance...37	
Commissioning	10
Connect MP3 player	25
Copy protection	13

D

Decommissioning.....	11
Definitions	5
Description	16
Device description.....	16
Device operating time	10
Dimensions	44
Directives	12
Disinfection	27
Disposal	11

E

Electrical installation	9
Electrical specifications.....	43
Environmental	11
Equipment.....	15

Error codes	27
Exhaust air ventilation.....	9
Export.....	12

F

Faults	27
Filter	
Clean	41
Cleaning intervals.....	29
Foreseeable misuse.....	7
Functions	23

G

Guarantee	14
-----------------	----

I

Information for the operator	7
Instruction signs	13
Intended use	7

L

Lamp sets	45
Lamps	
Control.....	45
Power levels	45
Rated power	45
Light sensitivity.....	19

M

Maintenance	26, 27, 30
Maintenance plan.....	31
Measurements	44
Medications.....	19
MP3 music	13

N

Noise level	43
-------------------	----

O

Operating altitude.....	10
Operating overview	22
Operator's obligations	7
Order numbers.....	45, 46

P

Packaging	12
Personnel training	8
Plastic surfaces (maintenance)	28
Play rights	13
Plug-in connector	9
Power output	43
Protective measures when performing maintenance work ..	26

R

Recycling	11
Removing the acrylic glass panel	33
Removing the canopy panel	34
Removing the front panel	36
Replacing the lamps	37
Replacing the starter	38
Risk of injury	19

S

Safety	9
Safety equipment	14
Scope of delivery	8, 15

Service kit	32
Setup	9
Setup location	10
Sickness	19
Skin sensitivity	20
Skin type	20
Spare parts Order numbers	46
Start application	24
Starter	45
Start-up settings	24
Storage	11
Symbols	6

T

Technical data	43
Time control	10
Transport	9
Treatment programmes	18

V

Ventilation	9
-------------------	---

W

Warranty	14
----------------	----

BEAUTY FLAIR

JK-SALES GMBH

Koehlershohner Strasse
53578 Windhagen/GERMANY
Phone: +49 (0) 2224/818-0
Fax: +49 (0) 2224/818-500
Internet: www.ergoline.de
E-Mail: info@ergoline.de

Ergoline